

Hydraulic ENT chair

Range 512 **51266T-ORL**







TECHNICAL DETAILS

- 3 sections ENT chair with hydraulic height adjustment
- Rotation at 360° (possibility to block the rotation)
- Synchronised backrest and legrest adjustment, with gas jack
- Adjustable head rest specialy developed for ENT exams
- Trendelenburg position
- Ergonomic upholstery with good density 30kg/cm3, M1 fire retardant class, thickness 8cm (specialy developed for an optimal comfort of the patient, with a lumbar support)
- Equipped with 2 removable armrests and a foot rest
- Chromed base

DIMENSIONS

- Width 55cm
- Adjustable height from 52 to 72cm
- Load capacity: 130kg in dynamic



Adjustable head rest specialy developed for ENT exams

ACCESSORIES




ACC514-ASSISE ACC514-JAMBIÈRE

Plastic protection for seatrest or legrest



ACC 993 4 retractable wheels, diam.50mm

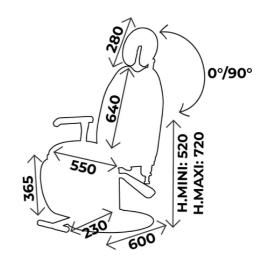





Warranty: 2 years (1 year for the upholstery)

UPHOLSTERY'S COLOURS

Fire retardant upholstery, M1 class

For the upholstery colours, you can choose between the Classic and Urban range Be careful! Urban colours with extra value

Upholstery: Classic range



Upholstery: Urban range



BASE COLOURS



MAINTENANCEGUIDE

DAILY MAINTENANCE

Upholstery: Cleaning with soapy water only (no solvents), can be disinfected with products such as Phagosept Spray, Salvanios pH10 at 0.75% and Anios Surfaces Hautes, without abrasive material

Base: cleaning with soapy water only (no solvents) without abrasive materials

WARRANTY

Mechanical/electrical part: 2 years

Upholstery: 1 year

CONTACT



Parc d'entreprises Visionis 278 rue de l'industrie 01090 GUÉREINS- FRANCE

COMMERCIAL SERVICE

Tel: +33(0)4 74 66 69 49

Fax: +33(0)474666950

export@carinamedical.com