



Picture of the model with all around controller (Ref. 99926)

Physiotherapy table

999 range

999 24

CARINA
Medical




CE
ISO 9001
ISO 13485

TECHNICAL DETAILS

- 2 sections physiotherapy table, with electric height adjustment
- Hand remote control
- Adjustment of the backrest in every position with a gas jack
- Stretched upholstery with good density 40g/cm², M2 fire retardant class, thickness 7 cm
- Height adjustment with an electric jack of 6000N capacity
- Face hole to give a better confort to patients, with plug cushion
- Equiped with exam sheet roller on the foot side
- Epoxy covered base, silver color in standard
- With 4 removable wheels, diam.35mm

DIMENSIONS

- Width 600mm
- Adjustable height from 48 to 95cm
- Load capacity: 150kg

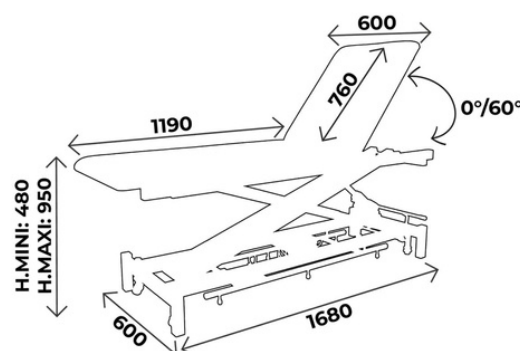

You may also like :

Réf. 99921 : Same as 99924 but 4 fixed feets

Réf. 99922 : Same as 99924 but with pedal control and 4 fixed feets

Réf. 99923 : Same as 99924 but with all around controller and 4 fixed feets

Réf. 99925 : Same as 99924 but with pedal control

Réf. 99926 : Same as 99924 but with all around controller



WHAT YOU NEED TO KNOW

- ✓ Delivery time: 15 days upon receipt of order
- ✓ Warranty: 2 years (1 year for the upholstery)

ACCESSORIES



ACC 538
Pair of clamps diam.14mm
width 600mm




ACC964
Pair of legrests on knee system
in black polyurethane,
diam.14mm, without clamp



ACC 76C
Unit clamp
to fix infusion stand diam.16mm




ACC 110
Blood test splint in colour
without clamp, the unit

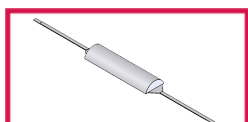

ACC 54001
Half moon headrest width 600mm
fire retardant M1 class



ACC 514-UNIVERSE
Plastic protection cover

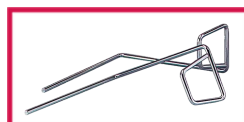

ACC 105
Pair of stainless steel stirrup
diam. 14mm, without clamp



SLFEMIREPOSE
Cover for legrests ACC964
in upholstery, the unit



ACC32012
Stainless steel infusion stand
diam. 16mm
2 plastic hooks, without clamp



ACC 110POLY
Polyurethane blood test splint
without clamp, the unit



ACCTETERGO
Ergonomic head rest

UPHOLSTERY'S COLOURS

Fire retardant
upholstery, M2 class

*For the upholstery colours, you can choose between the Ecomax, Classic and Urban range
Be careful ! Classic and Urban colours **with extra value***

Upholstery: ECOMAX range*



Galapagos



Bahamas



Blanc



Tahiti



Malouines

**Standard colors for ECOMAX, SIMPLEX and 999 ranges,
Classic and Urban colors available **on demand***

UPHOLSTERY'S COLOURS

Fire retardant
upholstery, M1 class

For the upholstery colours, you can choose between the Classic and Urban range

Upholstery: Classic range




Arno




Garonne



Potomac

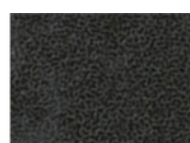


Amour

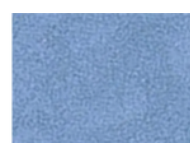

Mississippi



St Laurent

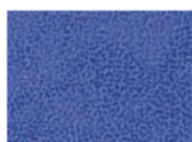

Rhône



Missouri



Seine




Amazon

Upholstery: Urban range



Abyss

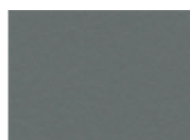

Cendre



Origan




Géranium




Régisse




Prune

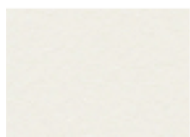


Perle

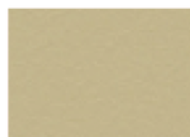

Ficelle

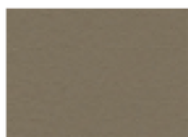


Terre

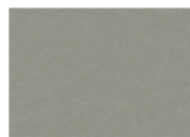


Béton



Chataigne

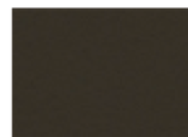

Chocolat

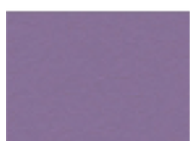

Pensée



Pétrole



Turquoise



Granny



Absinthe



Ananas



Orange



Fuschia



Tomate



MAINTENANCE GUIDE

DAILY MAINTENANCE

Upholstery: Cleaning with soapy water only (no solvents), can be disinfected with products such as Phagosept Spray, Salvanios pH10 at 0.75% and Anios Surfaces Hautes, without abrasive material

Base: cleaning with soapy water only (no solvents) without abrasive materials

WARRANTY

Mechanical/ electrical part: 2 years

Upholstery: 1 year



CONTACT



Parc d'entreprises Visionis
278 rue de l'industrie
01090 GUÉREINS- FRANCE

COMMERCIAL SERVICE

Tel: +33(0)4 74 66 69 49
Fax: +33(0)4 74 66 69 50
export@carinamedical.com